

Customer Details: Zimpli Medical Ltd  
Gladstone Street  
Greenbank Business Park  
Blackburn  
BB1 3ES

Customer Contact Name: Eejay Enyi  
Customer Email: eejay@zimplikids.com  
Purchase Order Number: Proforma

Date of Report: 22.04.20  
Date on Test: 17.04.20  
Date Received: 23.03.20

**MelBec Reference Number:** 16205  
**No. of samples:** 1

**Sample Details:**

Manufacturer/Supplier: Zimpli Medical Ltd  
Name of Product: Handsafe Alcohol Hand Sanitiser  
Storage: Room Temperature

**Experimental Conditions as outlined in the standard:**

Neutraliser: BU Broth  
Bacterial Strains: *E. coli* K12 NCTC 10538  
Product Dilution: Neat (as supplied)  
Test Product Application: 3ml applied to dry hands. Hands rubbed for 30s.  
Incubation Temperature: 36°C ± 1°C

**Conclusion:**

The test product was not inferior to the reference product and hence has met the requirements of EN 1500.

## Test Report for BS EN 1500:2013

**Results:**
**Inoculum (N)** – 4.90 x 10<sup>8</sup> cfu/ml

**Reference** Handrub: Reference (60% Propan-2-ol)

Subject	Hand	Mean Pre Value	Mean L+R	Mean Log Pre Value	Mean Post Value	Mean L+R	Mean Log Post Value	Log Reduction Factor
1	L	5.10E+06	4.55E+06	6.66	1.88E+02	2.89E+02	2.46	4.20
	R	4.00E+06			3.90E+02			
2	L	4.20E+05	6.80E+05	5.83	1.13E+03	9.00E+02	2.95	2.88
	R	9.40E+05			6.70E+02			
3	L	7.20E+06	8.70E+06	6.94	2.77E+03	2.52E+03	3.40	3.54
	R	1.02E+07			2.27E+03			
4	L	5.50E+06	4.18E+06	6.62	6.80E+02	2.29E+03	3.36	3.26
	R	2.86E+06			3.90E+03			
5	L	3.30E+06	3.90E+06	6.59	5.30E+02	5.75E+02	2.76	3.83
	R	4.50E+06			6.20E+02			
6	L	4.10E+06	4.30E+06	6.63	2.10E+02	3.15E+02	2.50	4.14
	R	4.50E+06			4.20E+02			
7	L	5.50E+06	7.05E+06	6.85	3.00E+03	2.58E+03	3.41	3.44
	R	8.60E+06			2.15E+03			
8	L	8.80E+06	8.75E+06	6.94	4.10E+03	4.30E+03	3.63	3.31
	R	8.70E+06			4.50E+03			
9	L	1.90E+06	1.01E+06	6.00	1.25E+02	1.93E+02	2.28	3.72
	R	1.20E+05			2.60E+02			
10	L	7.90E+05	7.80E+05	5.89	4.40E+02	8.35E+02	2.92	2.97
	R	7.70E+05			1.23E+03			
11	L	5.60E+05	7.10E+05	5.85	4.30E+02	1.35E+03	3.13	2.72
	R	8.60E+05			2.26E+03			
12	L	1.77E+05	1.78E+05	5.25	7.10E+01	1.75E+02	2.24	3.01
	R	1.79E+05			2.78E+02			
13	L	3.50E+05	4.25E+05	5.63	1.15E+03	1.31E+03	3.12	2.51
	R	5.00E+05			1.46E+03			
14	L	1.52E+06	8.20E+05	5.91	3.80E+02	7.30E+02	2.86	3.05
	R	1.19E+05			1.08E+03			
15	L	3.10E+06	3.80E+06	6.58	8.10E+02	1.16E+03	3.06	3.52
	R	4.50E+06			1.50E+03			
16	L	2.50E+05	2.95E+05	5.47	4.41E+02	4.26E+02	2.63	2.84
	R	3.40E+05			4.11E+02			
17	L	8.30E+05	1.17E+06	6.07	4.40E+02	3.45E+02	2.54	3.53
	R	1.50E+06			2.50E+02			
18	L	1.40E+05	3.05E+05	5.48	3.40E+01	1.97E+02	2.29	3.19
	R	4.70E+05			3.60E+02			

Imperial House, Kingsway, Haslingden, Rossendale, BB4 4QJ

info@melbecmicrobiology.co.uk Tel – 01706 214 492 www.melbecmicrobiology.co.uk

**Test Product: Handsafe Alcohol Hand Sanitiser**

Subject	Hand	Mean Pre Value	Mean L+R	Mean Log Pre Value	Mean Post Value	Mean L+R	Mean Log Post Value	Log Reduction Factor
1	L	9.20E+05	1.03E+06	6.01	1.99E+02	8.50E+02	2.93	3.08
	R	1.14E+06			1.50E+03			
2	L	7.40E+05	7.50E+05	5.88	7.32E+03	5.17E+03	3.71	2.16
	R	7.60E+05			3.02E+03			
3	L	6.20E+06	6.20E+06	6.79	6.30E+03	7.60E+03	3.88	2.91
	R	6.20E+06			8.90E+03			
4	L	2.30E+06	2.45E+06	6.39	2.00E+03	3.45E+03	3.54	2.85
	R	2.60E+06			4.90E+03			
5	L	2.40E+05	2.10E+05	5.32	9.50E+01	9.25E+01	1.97	3.36
	R	1.80E+05			9.00E+01			
6	L	1.64E+06	1.25E+06	6.10	1.13E+03	1.28E+03	3.11	2.99
	R	8.50E+05			1.42E+03			
7	L	3.70E+06	4.20E+06	6.62	3.00E+03	4.05E+03	3.61	3.02
	R	4.70E+06			5.10E+03			
8	L	4.80E+06	5.75E+06	6.76	1.39E+04	1.05E+04	4.02	2.74
	R	6.70E+06			7.10E+03			
9	L	5.50E+05	5.10E+05	5.71	1.54E+02	6.72E+02	2.83	2.88
	R	4.70E+05			1.19E+03			
10	L	1.58E+06	1.86E+06	6.27	4.40E+02	1.53E+03	3.18	3.09
	R	2.13E+06			2.61E+03			
11	L	1.23E+06	1.45E+06	6.16	1.17E+03	9.95E+02	3.00	3.16
	R	1.67E+06			8.20E+02			
12	L	7.83E+06	5.97E+06	6.78	8.50E+01	4.13E+02	2.62	4.16
	R	4.10E+06			7.40E+02			
13	L	1.31E+06	1.09E+06	6.04	1.17E+02	3.05E+02	2.48	3.55
	R	8.70E+05			4.92E+02			
14	L	6.30E+05	3.52E+06	6.55	2.95E+02	3.08E+02	2.49	4.06
	R	6.40E+06			3.20E+02			
15	L	3.80E+06	3.50E+06	6.54	2.57E+03	4.85E+03	3.69	2.86
	R	3.20E+06			7.12E+03			
16	L	5.00E+05	7.45E+05	5.87	7.40E+01	1.24E+02	2.09	3.78
	R	9.90E+05			1.74E+02			
17	L	2.20E+06	2.90E+06	6.46	8.70E+02	9.15E+02	2.96	3.50
	R	3.60E+06			9.60E+02			
18	L	1.74E+06	1.31E+06	6.12	2.20E+01	4.46E+02	2.65	3.47
	R	8.70E+05			8.70E+02			

Imperial House, Kingsway, Haslingden, Rossendale, BB4 4QJ

info@melbecmicrobiology.co.uk Tel – 01706 214 492 www.melbecmicrobiology.co.uk

**Summary:**

	<b>Mean of Pre-Values (Log)</b>	<b>Mean of Post-Values (Log)</b>
Reference Product	6.18	2.86
Test Product	6.24	3.04

	<b>Mean Log Reduction</b>
Reference Product	3.32
Test Product	3.20

	RP-PP	PP-RP
Mean lg R RP	3.59	3.04
Mean lg R PP	2.89	3.51
Difference of Means	0.70	-0.47

Absolute Difference of Differences (0.70 – (-0.47)) = 1.17  
 (<2.0 hence acceptable).

**Validations:**

Neutraliser Control B (NV<sub>B</sub> 3.0 x 10<sup>4</sup> – 1.6 x 10<sup>5</sup> cfu/ml) – Inoculum 9.80 x 10<sup>4</sup>cfu/ml

<b>Organism</b>	<b>Cfu/ml</b>
<i>E. coli</i> K12	7.6 x 10 <sup>2</sup>

Method Validation C (NV 3.0 x 10<sup>2</sup> – 1.6 x 10<sup>3</sup> cfu/ml) – Inoculum 3.27 x 10<sup>2</sup>cfu/ml

<b>Product</b>	<b><i>E. coli</i> (cfu/ml)</b>
Test Product 16205	3.30 x 10 <sup>1</sup>

**Verification:**

A complete set of results from at least 18 volunteers is available.

The overall means of the Ig prevalues for RP and PP is at least 5.00.

There are no more than three individual Ig reductions less than 3.00 occurring in RP.

The absolute difference of mean differences between Ig reductions of RP and PP of group RP → PP and group PP → RP is less than 2.00.

Control of weighted mean counts: quotient  $>5 < 15$ .

N is between  $1.5 \times 10^8$  and  $5.0 \times 10^8$

NV is between  $3.0 \times 10^2$  and  $1.6 \times 10^3$

NV<sub>B</sub> is between  $3.0 \times 10^4$  and  $1.6 \times 10^5$

B is  $\geq 0.0005 \times NV_B$

C is  $\geq 0.5 \times NV_0$

**Reduction Factor Difference:**

Subject	Reduction Factor Difference
1	1.11
2	0.72
3	0.63
4	0.41
5	0.48
6	1.15
7	0.42
8	0.57
9	0.84
10	-0.11
11	-0.44
12	-1.15
13	-1.04
14	-1.01
15	0.66
16	-0.94
17	0.03
18	-0.28


**Hodges-Lehmann 97.5% Upper Confidence Limits:**

Sorted Differences	Mean Pairwise Differences									
	1.15	1.11	0.84	0.72	0.66	0.63	0.57	0.48	0.42	0.41
1.15	1.15									
1.11	1.13	1.11								
0.84	0.99	0.98	0.84							
0.72	0.93	0.92	0.78	0.72						
0.66	0.90	0.89	0.75	0.69	0.66					
0.63	0.89	0.87	0.73	0.67	0.64	0.63				
0.57	0.86	0.84	0.70	0.64	0.61	0.60	0.57			
0.48	0.81	0.79	0.66	0.60	0.57	0.55	0.52	0.48		
0.42	0.78	0.77	0.63	0.57	0.54	0.52	0.50	0.45	0.42	
0.41	0.78	0.76	0.62	0.56	0.53	0.52	0.49	0.44	0.42	
0.03	0.59	0.57	0.43	0.37	0.34	0.33	0.30	0.25		
-0.11	0.52	0.50	0.36	0.30	0.27	0.26	0.23			
-0.28	0.43	0.42	0.28	0.22	0.19	0.18				
-0.44	0.35	0.34	0.20	0.14	0.11					
-0.94	0.10	0.09								
-1.01	0.07	0.05								
-1.04	0.05	0.04								
-1.15										

The median is between the 9th and 10th value:  $(0.42 + 0.41) / 2 = 0.42$

The mean pairwise differences that do not exceed the median (here 0.42) are computed. The critical values for Wilcoxon's matched pairs signed ranks test the entry for  $n=18$  and a one-sided 0.025 level of significance, the critical value of 40 is found. Hence,  $c = 40 + 1 = 41$ . The 41st entry is 0.57. Hence, the Hodges-Lehmann upper one sided 97.5% confidence limit for the difference in lg reductions between RP and PP is 0.57 which is less than the agreed inferiority margin of 0.6 lg units. Therefore, the hypothesis of inferiority of PP (test product) is rejected and it can be concluded that the test preparation PP is not inferior to RP (reference product).

Report Authorised By:

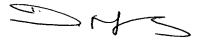

Dawn Mellors

-----  
Technical Director

This is the end of the test report